

SAPIENS-200: EPA 604 Phenols

Column: SAPIENS-200 P/N TR-580232


Dimensions: 30 m x 0,25 mm x 0,25µm

Injection: 1µl EPA 604 Phenols mix, 250 °C Split 1:50

Carrier gas: H₂, ct pressure 12 psi

Oven temperature: 50 °C (2 min) to 200 °C (5 min) @ 6°C/min

Detector: FID, 250 °C

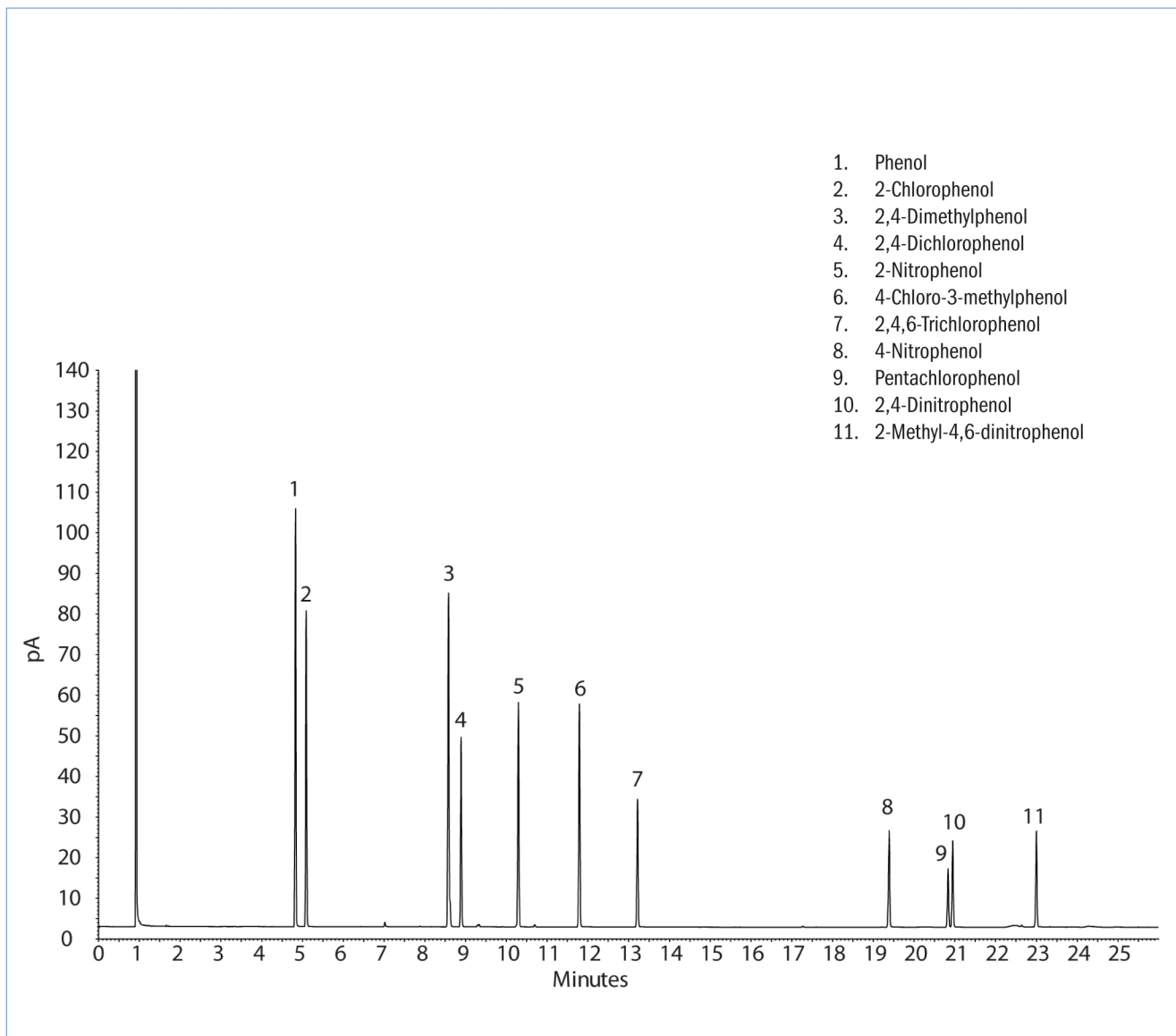

SAPIENS-200 TKG 1307

www.teknokroma.es

Teknokroma Analitica S.A.- Camí de Can Calders, 14 - 08173 Sant Cugat del Vallès (Barcelona) Spain
Tel: +34) 936 748 800 Fax: +34) 936 752 405 e-mail: comercial@teknokroma.es